



2018 Catalogue



 **major
science**

Innovative Life Sciences Tools

Chemical Engineering

Thermostirrer

The Major Science Thermostirrer is equipped with built-in stirrers on the bottom of each well. Combined with excellent temperature control, the instrument is well-suited for a variety of applications, including synthesis of chemical compounds, denaturation, derivatization, combinatorial chemistry, sample concentration, enzyme assays and process optimization.

Thermostirrer

▣ Thermostirrer

99-100

Thermostirrer

Remote Thermostirrer Module



*All images are for reference only, actual products might differ from the images.

*Please visit our website at www.majorsci.com for more product selection and detailed information.

*Technical specification subject to change without notice.

Thermostirrer Remote Thermostirrer Module

TS-8W
CM-M0512-001

Major Science Thermostirrers feature leading-edge speed performance and the broadest temperature control range to cover a variety of chemical and biological applications including synthesis of organic compounds, derivatization, drug solubility, protein crystallization and process optimization. Mixing efficiency is maintained uniformly across the sample wells due to independent magnets under each well and thus ensures consistency of experimental results. Furthermore, Major Science designed a remote thermostirrer module for studies that require a thermostatic mixer to function as a part of an automated process such as a polymer viscometers.

The Remote Thermostirrer Module allows you to set your own program with modbus PC connection and integrate with your own manufacture equipment. For more information with this module, please contact Major Science.




CE
TS-8W



User temperature calibration



Cooling System Connection Port

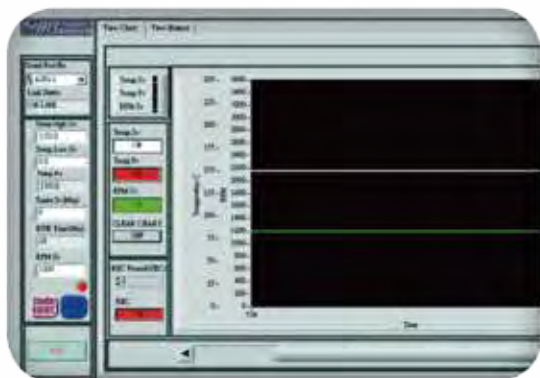

Functional control software package available



CE

CM-M0512-001
Remote Thermostirrer Module

Features:

- Microprocessor control with digital performance
- Brushless motor for individual well agitation
- Outstanding temp. control performance up to 200°C
- User temperature calibration
- Customization available
- Well-insulation around casework
- Over temperature protection
- LCD screen & timer (standalone model only)
- Highest stirring speed, 3500rpm in the market



Specifications

Cat. No.	TS-8W-110	TS-8W-220	CM-M0512-001
Controller	Digital microprocessor controller		
Display	2.6" LCD monochrome display		N/A
Motor	Brushless motor		
Number Of Position	8 wells (2x4) with individual magnetic stirrer or customized		6 wells with individual magnetic stirrer **
Well Diameter	Ø29.5mm, 62mm depth or customized		Ø23mm, 44mm depth **
Stirring Speed	500 - 3500rpm		500 - 2500rpm
Temperature Control Range	5°C above ambient to 200°C		
Temperature Increment	± 0.1°C		
Temperature Uniformity	± 0.7°C @ 150°C		N/A *
Temperature Accuracy	± 0.2°C @ 150°C		± 0.5°C @ 120°C
Temperature Calibration	Yes		
Operating Temperature	Ambient to 40°C		
Timer	99 (hr): 59 (min) / Continuous		
Block Material	PTFE coating with individual magnetic		Anodic treatment without PTFE coating **
Data Logging	RS 232 (Max. 2.5 meter long)		
Safety	Well-insulation around casework		
	Leakage proof for heating chamber		
	Over temperature protection		
Rated Voltage	110V~; 50/60Hz, 6.3A	220V~; 50/60Hz, 3.15A	110V~; 50/60Hz, 6.3A 220V~; 50/60Hz, 3.15A
Dimension (WxLxH)	Approx. 9.1"x11.8"x6.3" (230x300x160mm)		Approx. 5.2"x11.8"x6.3" (132.5x299.5x158.9mm)
Weight	Approx. 18.7lb (8.5kg)		Approx. 11lb(5kg)
Power	600W		
Vessel Diameter	Max. 28.5mm or customized		Max. 22mm **
Communication Port	RS-232		Modbus RTU by RS-232

* Temperature uniformity would be dependent on the integration to your automation system.

** Customized Specification.

Ordering Information

Cat. No.	Product Description
TS-8W-110	Thermostirrer 8 wells, 110V~
TS-8W-220	Thermostirrer 8 wells, 220V~
CM-M0512-001	Remote Thermostirrer Module 6wells, 220V~ / 110V~

Accessories

Cat. No.	Product Description
TS-8W-SW	Functional control and data logging software for Thermostirrer series. Including software CD and cable.



US Office

19959 Sea Gull Way
Saratoga, CA 95070
U.S.A.
T/ +1-408-366-9866
F/ +1-408-446-1107

Taiwan Office

No.19, Ln. 207, Huakang St.
Bade Dist., Taoyuan City 33464
Taiwan
T/ +886-3-3623319 / +886-2-2298-1055
F/ +886-3-3623133 / +886-2-2299-7871

Shanghai Office

Room 612, International business exhibition center,
9300 Hunan Road, Pudong, Shanghai, China
National toll-free No.:400-823-9177
T/ +86-21-50795277
F/ +86-21-50795277

India Office

D. No. 12-13-99, Satguru Apts. Extn.
Street. No. 3, Lane No.1, Tarnaka
Secunderabad – 500 017. India
T/ +91-40-27001515
T/ +91-40-27001586



www.majorsci.com
info@majorsci.com

